

Breeze® Self-Adhesive Resin Cement

Presenting Breeze® Self-Adhesive Resin Cement, one of the strongest*, most convenient, auto-mix self adhesive cement available.

When it's time to cement indirect restorations clinicians demand a product that is simple and convenient to use while providing high bond strength. After devoting time and energy creating a restoration to improve a patient's smile and self-esteem, don't settle for less than ideal adhesive materials for placement. **Make cementation a breeze.**



Breeze® Self-Adhesive Resin Cement

This self-etch formulation incorporates new resin technology for enhanced bond strength with greater reliability. Upon application, Breeze Cement quickly goes to work to condition dentin, enamel, and the restoration all in one quick and simple step. By eliminating the etching, priming and bonding steps, both postoperative sensitivity and procedure time are drastically reduced, improving peace of mind and freeing up precious chair time.

- > **Exceptional bond strength** - Adheres to dentin and enamel without a bonding agent
- > **Dual-cure** - Light-cure or self-cure for maximum versatility
- > **Convenient auto-mix delivery system** - Single-step accurate placement with minimal waste
- > **Contains fluoride** - Provides a sustained fluoride release to minimize the potential for caries
- > **Versatile shade selection** - Choose from A2, Translucent, and Opaque White to satisfy all cosmetic needs

1 > Dispense



Dispense Breeze® Cement directly into restoration.

2 > Place



Seat restoration.

3 > Cure




Light cure or self cure.

Photos courtesy of Dr. Michael Maroon

Achieve outstanding clinical results from start to finish in three easy steps.

Breeze Self-Adhesive Resin Cement is engineered to simplify the process of placing indirect restorations with reliable and consistent bond strength. Forget about etching, priming, bonding, hand mixing, or triturating capsules!

Bond Strength to Dentin†

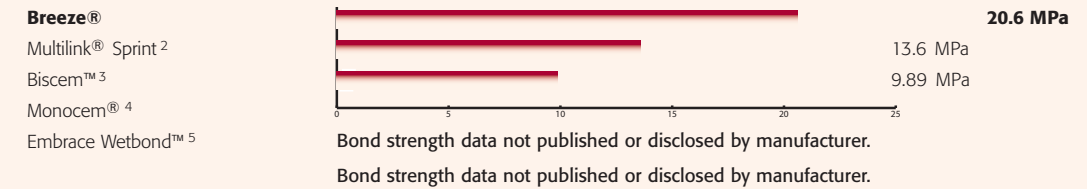

1 Manufacturer's published data 2 Multilink is a registered trademark of Ivoclar Vivadent, Inc. 3 Biscem is a trademark of Bisco, Inc. 4 Monocem is a registered trademark of Shofu Dental Corp. 5 Embrace Wetbond is a trademark of Pulpdent Corp.

What doctors are saying about Breeze Cement:

"Breeze® Self-Adhesive Resin Cement makes cementation easy whether it's a post or crown and clean-up is a snap. Whether you're luting metal, fiber, composite or ceramic restorations, Breeze® is my luting agent of choice"
- Gregori M. Kurtzman, DDS, MAGD, FACD, DICOI

"I'm a big fan of self etching resin cements and I'm especially a fan of Breeze® Cement. It is easy to use, easy to clean up, and provides an incredible bond to both metals and ceramics."
- John C. Flucke, DDS

"Breeze® Cement provides me with controllable working time and excellent bond strength. It has become my only choice for a self-adhesive resin cement."
- Michael Maroon, DMD
Founder & Fellow, Academy of Comprehensive Esthetics



One Cement – Multiple Applications

Breeze Self-Adhesive Resin Cement is formulated to accommodate virtually all cementation procedures.† Its advanced resin technology and low expansion ensure compatibility with multiple substrates including composite, semi-precious metals, precious metals, and cementable ceramics.



Photo courtesy of Dr. Michael Maroon

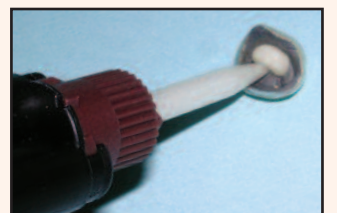


Photo courtesy of Dr. Gregori M. Kurtzman

| | Inlays/Onlays | Crowns | Bridges | Posts |
|--------------------------------------|---------------|--------|---------|-------|
| Metal/Porcelain Fused to Metal | X | X | X | X |
| Cementable Ceramic | X | X | X | X |
| Composite/Fiber Reinforced Composite | X | X | X | X |

*Data available upon request.

† Pentron Clinical recommends the Lute-It® Esthetic Luting System for cementing veneers and all-ceramic restorations that require bonding.

Breeze®

Self-Adhesive
Resin Cement

Ordering Information

Breeze® Self-Adhesive Cement Intro Kit REF N97

- 5 – 4ml/7.8gm syringes in the following shades:
A2 (3), Translucent, and Opaque White
- 1 – Bottle Silane 3ml/2.5gm
- 40 – Auto-mix Tips
- 10 – Intra-oral/Auto-mix Tips
- Instructions and Accessories

Breeze® Cement Refills

| | |
|--------------|----------|
| A2 | REF N97B |
| Translucent | REF N97A |
| Opaque White | REF N97E |

All cement refills include 1 - 4ml/7.8gm syringe of cement, auto-mix tips, and instructions.

Accessories

| | |
|--------------------------------|-----------|
| 1 – Bottle Silane 3ml/2.5gm | REF N04EA |
| 50 – Auto-mix Tips | REF N33CA |
| 100 – Auto-mix Tips | REF N33C |
| 50 – Intra-oral/Auto-mix Tips | REF N72AA |
| 100 – Intra-oral/Auto-mix Tips | REF N72AB |



For more information on our entire line of dental products visit us at www.pentron.com or call us at **1.800.551.0283**

Breeze®

Self-Adhesive
Resin Cement

- No etching
- No priming
- No bonding

No compromise.



800.551.0283 | www.pentron.com | 203.265.7397

Lute-It® | Bond-1® | Bond-It® | Nano-Bond® | Breeze® Adhesives Impressions Restoratives Crown & Bridge

Copyright © 2011 Pentron Clinical. All rights reserved.

BNLT_Rev06

